
FitScale W8 BT

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EN – USER MANUAL

The TrueLife FitScale W8 BT Wireless Smart Scale measures body proportions using bioelectrical impedance analysis (BIA). Technology and BIA are used to measure body fat and obtain total body composition data. The scale is intended for home use only.


Description of the scale:

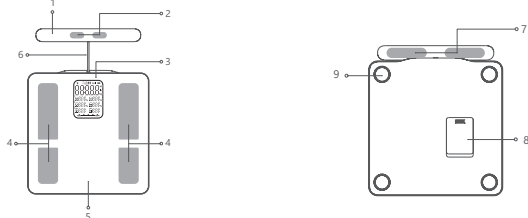1. Handle
2. Electrodes for thumbs
3. Display
4. Electrodes for feet
5. Body of the scale

6. Cable
7. Electrodes for fingers on the hand
8. Battery compartment cover
9. Rubber feet

Display Description:

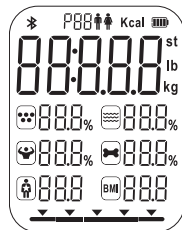	Proportion of adipose tissue		Hydration of the organism
	Proportion of muscle tissue		Proportion of bone tissue
	Body value		BMI

	Bluetooth icon		Battery level
	User code		Icon of a woman
	Icon of a man	lb	Weighing unit
kg	Weighing unit	kcal	Calorie unit
st lb	Weighing unit	%	Percentages


Body weight data:

	Thin
	Slim
	Normal weight
	Slight overweight
	Severe overweight

Inserting batteries:

Batteries are already inserted in the battery compartment. To use the scale, pull the protective foil out of the battery compartment.

If the low battery icon lights up, the batteries need to be replaced. Remove the battery compartment cover and insert the batteries into the compartment. When inserting, observe the correct polarity, which is marked on the bottom of the compartment. Put the cover back on.

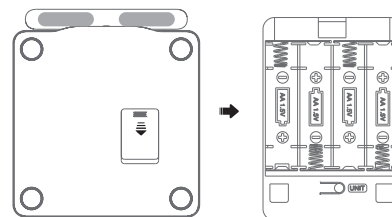Caution: Dispose of discharged batteries in an environmentally friendly manner. Do not combine different types of batteries, old and new batteries, etc.

Note: If you will not be using the scale for a long period of time, remove the batteries.

Note: After inserting the batteries, or if you move the scale, the weight will be displayed. Wait approximately 6 seconds for the scale to switch off. Then step lightly on the scale to turn it on. When "0.00 kg" is displayed, the scale is ready for use.

App installation:

Connect the FitScale to your mobile phone or tablet to measure all body proportions, such as subcutaneous fat. Connect the TrueLife FitScale W8 BT to your mobile phone via Bluetooth. To do this, you need a mobile app. The easiest way to get

the mobile app is to download it from the Appstore/Google Play or by scanning the QR code of the box or manual. Follow these steps for installation:

<https://www.trueLife.eu/01570/app>

Scan the QR code or search for the "MovingLife" app in the Apple Store/Google Play and download.

Sign up

2. Pairing the scale with the app



Turn on the Bluetooth function on your device. Stand on the scale. The scale automatically switches on when you step on its surface.


Pair the app with the scale

Weighing:

1. Place the scale on a firm and level surface. Open the app.
2. Press the scale to turn it on. Stand with your bare feet on the electrodes on the left and right side of the scale. Grasp the handle with both hands. The thumbs of the hands should rest on the electrodes on the top of the handle, the fingers of the hands on the electrodes on the bottom of the handle. Straighten up and bend forward so that your arms are at an angle of approximately 45° with your body.

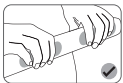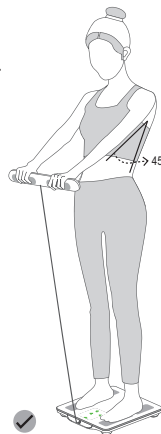3. When the weighing value on the display stops blinking and lights up, the analysis of the body values starts. Once completed, the complete analysis is displayed on the screen and transferred to the apps on your phone. The data is stored here.

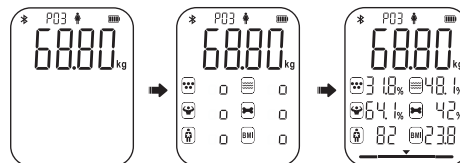Tips:

The entire weighing and analysis process takes approximately 20 seconds.


If the scale is not added to the app before weighing, only the body weight will be shown on the display, the analysis will not start. We recommend adding the scale to the app before weighing.

Note: Due to the continuous development of the application, the illustrations in the manual may differ slightly from those in the application. Final and accurate information is displayed in the app.

Note: The scale supports adding multiple users. A maximum of 15 groups of historical data can be stored for each user. When the app is turned on, the data will be automatically synced to your smartphone.

Tips for proper weighing:

- For consistent and accurate measurement results, please use at the same time of day under approximately the same conditions if possible. Avoid weighing in these cases:
 1. During or immediately after strenuous exercise
 2. Right after a shower, bath or sauna
 3. After heavy alcohol consumption
 4. During or immediately after a meal
 5. During fever
- The weighing unit can be changed from kilograms (kg) to pounds (lb.) in the mobile app. The height data (meter/foot) will also be changed.
- Place the weight correctly on a flat surface or on a solid base. Reset the scale before the first weighing. The scale should show 0.00 kg on the display before use.

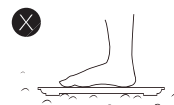Scale in profile on a flat surface and a scale on an uneven surface

- Weigh body fat and other body parameters correctly: stand on the scale upright, with both feet equally. The scale automatically switches on when you step on its surface. Be barefoot and stand upright. Weighing with socks will produce inaccurate results. Dirty feet can also affect weighing results.

Possible complications:

0	Measurements	Lo	Low battery
---	Error during weighing	Err	Overload
ErrL	Too low a body fat value	ErrH	Too high a body fat value

Part of the package:

Included with your TrueLife FitScale W8 BT scale is:

1. TrueLife FitScale W8 BT
2. 4x 1.5V AA batteries
3. User manual + safety regulations

Basic specifications:

Dimensions: 30 × 33.4 × 3.4 cm

Weight: 2,08 kg

Batteries: 4x 1.5 V AA

Weighing range: up to 180 kg

Weighing deviation: 0.1 kg

Units: kg/h lb/lb

PRODUCER RESPONSIBILITY

TrueLife products are warranted against manufacturing defects. In the event of such defects, elem6 s.r.o. is obliged to repair or replace the defective part or product with a comparable part or product at its discretion. Transfer, alteration or tampering with the product will exceed the scope of this warranty. This is the case even if the product is lost or damaged due to a defect, negligence or other error. The warranty does not cover damage resulting from use, accident or normal wear and tear. Elem6 s.r.o. assumes no liability for any accident, injury, death, loss or other claim related to or arising out of the use of this product. In no event shall elem6 s.r.o. be liable for incidental or consequential damages related to or resulting from the use of this product or any of its components. Returns or exchanges of replacement parts and/or products may be subject to shipping, handling, replacement and/or restocking charges. This product should be stored out of the reach of children. Use by children is only allowed under adult supervision.

EU Declaration of Conformity

elem6 s.r.o. hereby declares that the TrueLife FitScaleW8 BT complies with the essential requirements and other relevant provisions of Directive 2014/53/EU and 2014/30/EU. All TrueLife products are intended for sale without restriction in Germany, the Czech Republic, Slovakia, Poland, Hungary and other EU member states. The full Declaration of Conformity can be downloaded from www.truelife.eu/support/doc
Frequency bands in which the radio equipment operates: 2.4 ~ 2.48 Ghz

SAFETY WHEN WORKING WITH THE EQUIPMENT

1. Please read the enclosed safety instructions carefully.
2. For proper handling of the product, please read the current version of the manual, which is available for download at: www.truelife.eu.
3. Printing errors and information changes are subject to change.

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