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---

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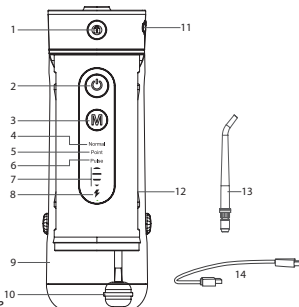
## PRODUCT DESCRIPTION

1. Button to release nozzle head
2. On/Off button
3. Mode change button
4. Normal mode icon
5. Point mode icon
6. Pulse mode icon
7. Water pressure diodes
8. Charging light
9. Water tank
10. Suction tube
11. Charging socket
12. Water tank cap
13. Interchangeable nozzle
14. Charging cable

**Warning:** There may be residual water in the device. This is distilled water used during quality control testing.

## INSTRUCTIONS FOR USE

Before the first use, remove all packaging material.

If the water flosser has been stored in a low temperature environment during transport, leave it at room temperature for at least two hours before use.

1. Connect the head to the shower. A click indicates correct insertion. You can choose from the following headers:
  - **Standard Head** – For daily plaque removal from interdental spaces and along gum lines
  - **Periodontal Head** – For cleaning periodontal pockets using a rubber tip
  - **Orthodontic Head** – For removing plaque from brackets and wires of braces and other hard-to-reach places
  - **Bubble Spraying Head** – For people with sensitive gums and people with braces
  - **Tongue Cleaner Head** – For effective cleaning of dirt on the tongue
2. Slide out the water tank. Open the water cap and fill the tank with water whose

temperature does not exceed 40°C. The tank capacity is 180 ml. Then close the lid and make sure that the water does not overflow. It is also possible to fill the water tank through the opening located in the rear upper part or by sliding the entire tank out of the body of the device.

3. Hold down the mode change button and select the desired shower mode. Press the same button to set the pressure level (1-4 levels). You can change the mode or pressure even while showering. Modes:

- **Normal mode** – suitable for daily cleaning. Pressure Range: 120/102/83/61 PSI
- **Point mode** – suitable for effectively removing food residues between the teeth, thanks to a program of three quick and precise sprays under high pressure followed by a short pause. Pressure Range: 67/56/46/37 PSI
- **Pulse mode** – suitable for massaging teeth and gums. Pressure Range: 120/114/100/83 PSI

4. Leave your mouth slightly open to allow water to rinse from your mouth into the sink. Press the On / Off button to start the flosser. Direct the water jet at a 90° angle to the teeth and gums. Slowly run the stream of water around your teeth. Direct the stream of water onto the gums and into the space between the teeth. If you have braces, you can direct the water between the joints of the braces.

5. When you have finished flossing, switch off the device using the On / Off button. Open the water tank cap and pour out any remaining water. Switch the device on again with the On / Off button to drain the water from the pump. Press the On / Off button to turn off the device and close the water tank cap.

**Warning:** To prevent the build-up of bacteria, ensure that no water remains in the device. If you do not intend to use the device for a longer time, make sure that it is thoroughly dried before storing it.

## CHARGING

The water flosser is not supplied with a fully charged battery. Before using for the first time, fully charge the water flosser.

Plug one end of the USB cable into the charging socket at the rear of the water flosser. Plug the other end of the USB cable into an adapter (adapter not included). Insert the adapter into the plug socket.

It takes approximately 4 hours to fully charge the water flosser.

**Warning:** To ensure maximum battery life, we recommend charging the water flosser once a week.

**Warning:** Charge the oral irrigator in a horizontal position on a stable surface. Keep the charging port dry and clean.

**Warning:** If you are not going to use the oral shower for a long time, we still recommend charging it regularly at least once every three months.

## STATUS OF CHARGE DETECTION

While the oral irrigator is off, press and hold the on/off and mode change buttons to check the battery charge status. The number of illuminated water pressure LEDs indicates the state of charge (4 LEDs indicate 75-100% charge; 3 LEDs indicate 50-75% charge; 2 LEDs indicate 25-50% charge; 1 LED indicates less than 25%).

**Warning:** When the charging light turns red for about 15 seconds, the charge level has dropped below 25% and the battery needs to be charged.

## INSERTING AND REPLACE THE NOZZLE HEAD

Insert the nozzle head into the hole at the top of the water flosser. You will hear a click when the nozzle head has been correctly attached.

To remove the nozzle head from the main body, press and hold the release button and then detach the nozzle head from the water flosser.

The nozzle head rotates 360 °.

**Warning:** Only insert and remove the nozzle head when the water flosser is switched off.

**Warning:** Nozzle heads are items of personal hygiene and are not covered by warranty.

**Tip:** You can store the head directly in the body of the oral shower. After removing the head, insert it into the hole located in the upper back part.

## MAINTENANCE

Use only water and neutral detergents to clean the device. Never use caustic or abrasive cleaners as they may damage the device. Do not use water whose temperature exceeds 50 °C. Keep the device away from high temperatures and avoid direct sunlight.

**Warning:** Do not submerge the device in water.

### 1. The device

Remove dirt by wiping with a towel or cloth with liquid soap and water, then ensure all water is rinsed off and wipe dry.

Do not immerse the device in water.

**Warning:** Water must not enter the USB socket.

### 2. Water tank

Wash the tank with water and wipe with a dry cloth.

Dry the tank thoroughly if you do not intend to use it for a longer time.

**Warning:** For greater ease of cleaning, the water tank can be removed from the body of the water flosser.

### 3. Nozzle head

Wash the nozzle head with water and wipe with a cloth.

### 4. Suction tube

Wipe the suction tube with a dry cloth.

Do not bend the suction tube, otherwise it may be damaged and the water flow interrupted during use.

## POSSIBLE COMPLICATIONS

| Problem                   | Cause                                         | Solution                                                      |
|---------------------------|-----------------------------------------------|---------------------------------------------------------------|
| The head cracks           | The head has been used for a long time        | Replace the head with a new one                               |
| The water tank is leaking | The seal on the body of the device is damaged | Contact the manufacturer with a request to provide a new seal |

|                                                        |                                        |                                                                |
|--------------------------------------------------------|----------------------------------------|----------------------------------------------------------------|
| Water is leaking from the connection point of the head | The head is not attached correctly     | Remove the head and reattach it when you hear a click          |
| Water pressure is very low                             | Oral shower is used upside down        | When using, the head must be placed higher than the water tank |
|                                                        | The water tank is empty                | Fill the tank with water                                       |
|                                                        | The water pressure level is set to low | Increase the water pressure level                              |

## SPECIFICATIONS

|                    |                                                                                |
|--------------------|--------------------------------------------------------------------------------|
| Battery life       | 12 days (2x daily 2 minutes)                                                   |
| Modes              | 3 (Normal, Point, Pulse)                                                       |
| Pressure levels    | 4 (for each mode)                                                              |
| Timer              | Yes (2 minutes)                                                                |
| Memory             | Yes                                                                            |
| Additional heads   | M-series (Standard, Orthodontic, Periodontal, Bubble Spraying, Tongue Cleaner) |
| Tank volume        | 180 ml                                                                         |
| Charging           | USB cable                                                                      |
| Pressure range     | 30-130 PSI                                                                     |
| Energy consumption | 10W                                                                            |
| Voltage            | 3.7V                                                                           |
| Battery capacity   | 1100mAh                                                                        |
| Water resistance   | IPX7                                                                           |
| Weight             | 236 g                                                                          |
| Dimensions         | 63×52×145 mm                                                                   |

## PRODUCT LIABILITY

TrueLife products are guaranteed to be free from manufacturing defects. In the event of such defects, elem6 s.r.o. will repair or replace the defective part or product with a comparable part or product at its sole discretion. The warranty will be deemed void if the device is altered or tampered with. The warranty does not cover damage resulting from use,

accident or normal wear and tear. elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss, or other claims related to or resulting from the use of this product. elem6 s.r.o. shall in no event be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Return or exchange of spare parts and / or products may be subject to shipping, handling, and replacement or replenishment fees.

This product should be stored out of the reach of children. Use by children should take place only under the direct supervision of an adult.

## **SAFETY GUIDELINES WHEN OPERATING THE DEVICE**

1. Read the enclosed safety instructions carefully.
2. The latest version of the manual, including instructions for operating the device, can be downloaded at: [www.truelife.eu](http://www.truelife.eu).
3. Printing errors and changes of information reserved.

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