# Care Q10 BT

Infrared thermometer

## CONTENTS

1. BOX CONTENTS.....	2
2. PRODUCT DESCRIPTION.....	2
3. INSTALLING THE BATTERIES.....	3
4. SETTINGS.....	4
5. GUIDELINES FOR MEASURING TEMPERATURE.....	4
6. TAKING A READING.....	5
7. DISPLAYING THE STORED READINGS.....	6
8. TRANSFERRING READINGS TO A MOBILE DEVICE.....	6
9. ERROR MESSAGES AND TROUBLESHOOTING.....	7
10. SPECIFICATIONS.....	7


## 1. BOX CONTENTS

1 × TrueLife Care Q10 BT thermometer	1 × Manual
2 × AAA alkaline battery	1 × Safety instructions
1 × Case	

The TrueLife Care Q10 BT uses infrared technology to achieve highly-accurate results. The Care Q10 BT thermometer measures body temperature in two possible ways, from the body surface, or with a special temperature sensor in the ear for even more accurate readings.


## 2. PRODUCT DESCRIPTION

- |                             |                        |
|-----------------------------|------------------------|
| 1. Display                  | 6. Sensor              |
| 2. Fever indicator          | 7. Change unit °C/°F   |
| 3. Power/Measurement button | 8. Sensor cover        |
| 4. Adult/Child mode button  | 9. Memory/Sound button |
| 5. Measurement mode button  | 10. Battery cover      |


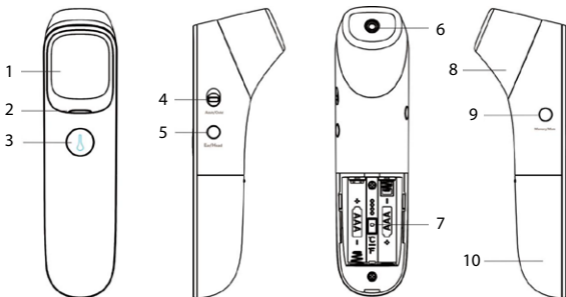
## Display description

1. Surface mode
2. Body mode
3. Child mode
4. Ear temperature mode
5. Mute
6. Measurement unit
7. Battery status indicator
8. Memory mode
9. Reading
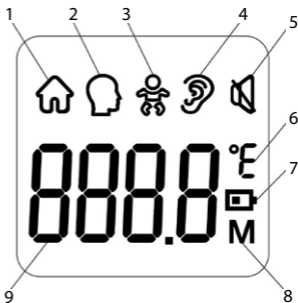
## 3. INSTALLING THE BATTERIES

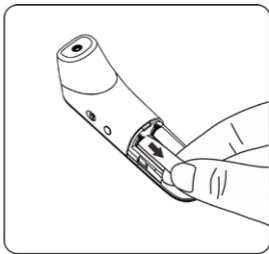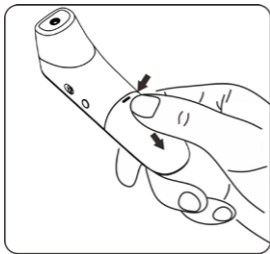Slide off the battery cover at the back of the thermometer and insert 2 x AAA batteries. Make sure that the batteries are inserted with the correct polarity as indicated in the battery compartment. Do not use batteries after their expiration date.

If the low battery indicator appears on the display after the thermometer is switched on, the batteries must be replaced. When battery power is critically low, the device will not take a reading.

If you will not be using the device for an extended period of time, remove the batteries.

## 4. SETTINGS

### Setting the temperature units

The thermometer reading can be displayed in degrees Celsius (°C) or Fahrenheit (°F). When the thermometer is turned off, remove the battery cover and select the unit of measurement using the switch between the batteries. Set the switch the up position for Celsius (°C) or to the down position for Fahrenheit (°F).

### Sound settings

The reading can be accompanied by an audible signal. Follow the instructions below to turn the beep on or off.

1. Switch on the thermometer with the power/measurement button.
2. Press and hold the memory/sound button for about 2 seconds.
3. When turning on the sound, the thermometer will beep. When turning off the sound, the mute symbol will appear on the display.

### Setting adult or child mode

The thermometer can measure temperature in adult or child mode, on the forehead or in the ear. Child mode is specially designed for children under 12 years of age.

Select the desired mode using the adult/child mode button. When the child mode is selected, the child symbol is shown on the display.

## 5. GUIDELINES FOR MEASURING TEMPERATURE

To ensure that the reading always shows your body temperature accurately, you must take into account the following factors that may affect the measurement.

1. It is important to know the normal healthy temperature of each individual. This is the only way to accurately diagnose a fever. To determine normal temperature, take several temperature measurements when the individual is healthy. Measure the temperature with a standard digital thermometer to confirm.
2. Users must be indoors for 30 minutes before taking a reading. Note: Users and the thermometer should be in the same temperature environment for at least 10 minutes before taking a reading.
3. Users should not drink, eat or be physically active before taking a reading, including bathing, showering, shampooing or drying their hair. Before each measurement, remove any head covering, remove hair from the forehead and wait 10 minutes.
4. Oils or cosmetics on the forehead may lead to a lower reading than the actual body temperature.
5. Keeping your hand on your forehead for any length of time will effect the reading.
6. Do not measure the temperature on scarred skin, open wounds or abrasions.
7. Do not measure the temperature on a sweaty forehead. The resulting reading may be distorted.

8. Do not take readings while breastfeeding or immediately after breastfeeding.
9. Do not use this thermometer outdoors.
10. Do not leave the thermometer near places that are very hot, such as fireplaces and stoves.
11. The thermometer sensor is the most sensitive part of the device. Do not touch the sensor. Reading accuracy may be affected if the sensor is damaged or dirty.
12. If the thermometer is stored in an environment with a significantly different temperature than where you will take the reading, keep it in the reading location for approximately 30 minutes before use.
13. It is not intended for use in an oxygen-rich environment or in the presence of a flammable anesthetic mixtures of air, oxygen or nitrogen dioxide.

## 6. TAKING A READING

The thermometer offers three temperature measurement modes.

- a) Forehead mode
- b) Ear mode
- c) Surface mode

Switch between the modes using the measurement mode button. You can identify the selected mode by the symbol on the display.

Taking a reading is simple. Follow the instructions below.

1. Press the power/measurement button to turn the thermometer on.
2. The currently set measurement modes are displayed and the thermometer is ready.
3. a) In forehead mode, aim the thermometer at the centre of the forehead at a distance of 1 to 3 cm and then press the power/measurement button.  
b) In ear mode, first remove the sensor cover. Then insert the sensor into the ear canal and press the power/measurement button.

**Caution:** Do not force the thermometer into the ear canal. The ear canal may get injured. Be careful when measuring the temperature of a child with a narrow ear canal.

- c) In surface mode, point the thermometer at the centre of the object you want to measure at a distance of 1 to 2 cm and press the power/measurement button.
4. The reading will be displayed.
5. The thermometer will automatically turn off if it is idle for more than 10 seconds to extend battery life.

### Colour and sound fever alert

The thermometer reading is accompanied by a colour and a sound when measuring on the forehead or in the ear.

When the temperature on the forehead or in the ear is between 35.0 °C/95.0°F and 37.5 °C/99.5°F, a long beep sounds and the indicator glows green, indicating that the body temperature is normal.

When the forehead temperature is between 37.6°C/99.6°F and 42.2°C/108.0°F and the in-ear temperature is between 37.6°C/99.6°F and 42.0°C/107.6°F, 3 short double beeps sound and the indicator glows red. The body temperature is high. The person may have a fever.

## 7. DISPLAYING THE STORED READINGS

The thermometer can store 20 readings. When the memory capacity is exceeded, the oldest stored entry will be overwritten.

Follow the instructions below to display the measured values.

1. Turn on the thermometer.
2. Press the memory/sound button. The memory mode symbol - the letter „M“ appears on the display.
3. Switch between the stored readings using the memory/sound button. Number 1 shows the last measurement, while number 20 shows the oldest.

## 8. TRANSFERRING READINGS TO A MOBILE DEVICE

The TrueLife Care Q10 BT thermometer allows you to transfer your readings to your mobile device so you can carry them with you at all times. The software automatically saves the measured data and draws a temperature curve graph to show the status of temperature changes in real time.

Below you will find instructions on how to transfer the data.

First, you need to download the app, install it and register.

1. Download the free app to your mobile device, which can be found at the link below. The app is designed for mobile devices or other devices running Android or iOS.
2. Install and open the application.
3. Proceed by registering in the app. Fill in your email address and password. This information serves as your login details.

[truelife.eu/56915/app](https://truelife.eu/56915/app)

## 9. ERROR MESSAGES AND TROUBLESHOOTING

If any of the following problems occur when taking a reading, please refer to the table below.

Problem	Probable Cause	Remedy
<b>Nothing appears on the display when the thermometer is turned on.</b>	The batteries are low.	Replace all batteries with new ones.
	The batteries are inserted incorrectly.	Check and, if necessary, insert the batteries correctly so that their polarity matches the markings in the battery compartment.
<b>Error message 1 (E1)</b>	The ambient temperature is not between 10 and 40°C (50°F ~ 104°F).	Place the thermometer in a room for at least 30 minutes at a room temperature between 10 and 40°C (50°F ~ 104°F).
<b>"Hi" message</b>	In forehead measurement mode: the temperature is higher than 42.0°C (109.4°F). In ear measurement mode: the temperature is higher than 42.2°C (109.4°F). In surface measurement mode: the temperature is higher than 100°C (212°F).	Carefully read the instructions on how to take a reading and then take another reading.
<b>"Lo" message</b>	In forehead measurement mode: the temperature is lower than 35.0°C (93.2°F). In ear measurement mode: the temperature is lower than 35.0°C (93.2°F). In surface measurement mode: the temperature is lower than 0°C (32°F).	Carefully read the instructions on how to take a reading, make sure the sensor is clean, and then take another reading.

## 10. SPECIFICATIONS

<b>Measuring range</b>	Body temperature: 35.0°C–42.2°C (93.2°F–109.4°F) Surface temperature: 0°C–100°C (32°F–212°F)
<b>Measurement deviation</b>	Body temperature: ± 0.2°C (0.4°F) Surface temperature: ± 1.0°C (2.0°F)
<b>Measurement distance</b>	Forehead body temperature: 2–3 cm Surface temperature: 1–2 cm
<b>Measurement time</b>	≤2 s

Memory	20 measurements
Battery	2 x 1.5V AAA
Operating conditions	Temperature: 10°C–40°C (50°F–104°F), Humidity: ≤ 95% RH
Storage conditions	Temperature: -25°C–55°C (-13°F–131°F), Humidity: ≤ 95% RH
Dimensions	155.2 × 39.6 × 49.1 mm
Weight	101 g (with batteries)

## SAFETY WHEN WORKING WITH THE DEVICE:

1. Read the attached safety instructions carefully.
2. For proper handling of the product, read the latest version of the manual, which is available for download at: [www.truelife.eu](http://www.truelife.eu).

Printing errors and information changes reserved.

*Note: This product is intended for use in households and similar environments. We declare that if the TrueLife Care Q10 BT is used as a professional medical device, a safety inspection (BTK) must be performed regularly at the following intervals: 24 months after the purchase of the product, then every 12 months.*

### Importer:

elem6 s.r.o.,  
Braškovská 15, 16100 Prague 6,  
Czech Republic

Made in P.R.C.

CE 0482 RoHS

Model: JPD-FR409-BT

MedPath GmbH  
Mies-van-der-Rohe-Strasse 8, 80807 Munich, Germany

Shenzhen Jumper Medical Equipment Co., Ltd.  
Address: D Building, No. 71, Xintian Road, Fuyong Street, Baoan, Shenzhen,  
Guangdong, China, 518103

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